


SELECTAIR®

Low Air Loss Mattress Replacement

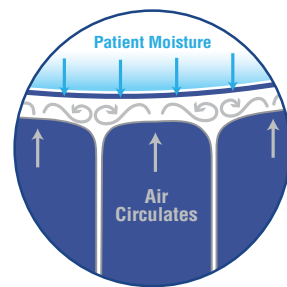

User-Friendly Control Unit

Rapid inflation, easy set-up and a motor so quiet, the patient will hardly know it's there!

The SelectAir® System

Quietly circulates air through individual cells and allows air to escape from the mattress.

Integrated Side Bolsters

Keeps patient centered in the bed and reduces the risk of entrapment.

Wicking Action:

Air escapes through tiny holes on the top of the mattress and circulates under a vapor-permeable cover to wick moisture away from patient.

SELECTAIR® Features

Firm Mattress Mode: Helps facilitate nursing care during dressing changed and patient repositioning.

Upright Patient Mode: Increases pressure to add support when the head of the bed is elevated.

Foot Pillow: Prevents unnecessary skin breakdown due to patient migration (optional).

CPR/Emergency Deflate: When engaged, allows quick deflation of the system in case of emergency.

Zoned, "true" low air loss used in acute care hospitals nationwide

Weight Limit: 400 lbs. (181 kg)

Mattress Size: Fits most standard bed frames. Custom sizes available. Inquire for details.

Cover: Water resistant, vapor-permeable, flame resistant and biocompatible

Warranty: One Year Limited Warranty on control unit and soft goods

Medicare Code: E0277 Group 2 Support Surface

LIT-0002 | R100121

 Manufactured in the USA by
MOXI
ENTERPRISES
Innovative health solutions

1939 S. Vandeventer Ave.
St. Louis, MO 63110
T 1-800-942-7678
E info@moxiusa.com
W moxiusa.com